

LIFEPAK CR° Plus defibrillator LIFEPAK EXPRESS° defibrillator

OPERATING INSTRUCTIONS



LIFEPAK CR®Plus defibrillator LIFEPAK EXPRESS® DEFIBRILLATOR

OPERATING INSTRUCTIONS

Responsibility for Information

It is the responsibility of our customers to ensure that the appropriate person(s) who may use this defibrillator have access to the information in this manual, including general safety information provided in Section 1.

Device Tracking

The U.S. Food and Drug Administration requires defibrillator manufacturers and distributors to track the location of their defibrillators. If the device is located somewhere other than the shipping address or the device has been sold, donated, lost, stolen, exported, destroyed, permanently retired from use, or if the device was not obtained directly from Physio-Control, please do one of the following: register the device at http://www.physio-control.com, call the device tracking coordinator at 1.800.426.4448, or use one of the postage-paid address change cards located in the back of this manual to update this vital tracking information.

Revision History

These operating instructions describe the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators with software Version 3.0 or later.



LIFEPAK, LIFEPAK CR, LIFEPAK EXPRESS, and QUIK-COMBO are registered trademarks of Physio-Control, Inc. ADAPTIV, CODE-STAT, CHARGE-PAK, QUIK-PAK, and Shock Advisory System are trademarks of Physio-Control, Inc. IrDA is a registered trademark of Infrared Data Association. Specifications are subject to change without notice.

©2009 Physio-Control, Inc. All rights reserved

Publication Date: 1/2009



TABLE OF CONTENTS

1 Introduction	
About Automatic External Defibrillators	1-2
Indications for Use	1-2
Contraindications	1-2
Why the Need for Defibrillators	
Terminology	
Text Conventions	
Safety Information	
Safety Terms	
General Warnings and Cautions	
Symbols	
About the LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillators	1-9
Capabilities and Features	
2 Getting Started	
Unpacking and Inspecting Your LIFEPAK CR Plus or LIFEPAK EXPRESS Defib	rillator2-2
Where to Locate Your LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator	
Controls, Indicators, and Labels	
Outside Controls, Indicators, and Labels	
Inside Features	
3 Using the Defibrillator	
Warnings and Cautions	
Responding to a Sudden Cardiac Arrest emergency	
Basic Steps for Using the LIFEPAK CR Plus or LIFEPAK EXPRESS Defibi	
What to Do After Emergency Medical Personnel Arrive	
What to Do After Using Your Defibrillator	
Voice Instructions and Tones	
Troubleshooting	

4 Data Storage	
Overview of Data Storage	4-2
Data Stored by Your Defibrillator	
Test and Service Data	4-2
Event and Test Log	4-3
5 Caring for Your Defibrillator	
Maintaining a State of Readiness	
Storing Your Defibrillator	
Cleaning Your Defibrillator	
Replacing the CHARGE-PAK Battery Charger and the QUIK-PAK Electrode Packet	
Replacing the CHARGE-PAK Battery Charger	
Replacing the QUIK-PAK Electrode Packet	
Obtaining Authorized Service	
Recycling Information	
Recycling Assistance	
Preparation	
Recycling of Disposable Electrodes	
Packaging	
Supplies, Accessories, and Training Tools	
<u> </u>	5-0
6 Defibrillator Operating Settings	
Operating Settings and Setup Configuration	6-2
A Specifications	
B Shock Advisory System	
C User's Checklist	
D Electromagnetic Compatibility Guidance	

Index

INTRODUCTION

This section provides background information about defibrillation and includes an overview of LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillator features.

About Automatic External Defibrillators	page 1-2
Safety Information	1-4
Symbols	1-7
About the LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillators	1-9

ABOUT AUTOMATIC EXTERNAL DEFIBRILLATORS

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are automated external defibrillators (AEDs). For many years, defibrillators have been used only by medical professionals to treat victims in sudden cardiac arrest (SCA). Today, the ability of defibrillators to save lives is so widely recognized that people once trained to do only cardiopulmonary resuscitation (CPR) can now use defibrillators.

After electrode pads are applied to the victim's chest, the defibrillator analyzes the victim's heart rhythm. If a shockable rhythm is detected, the defibrillator will either deliver an intense pulse of electricity (shock) to the heart muscle (fully automatic model) or direct the responder to deliver the shock (semiautomatic model). The defibrillator delivers shocks through the electrode pads on the victim's chest.

When this pulse of electricity is delivered, it is called defibrillation. Defibrillation is recognized for treating life-threatening heart beat irregularities, such as ventricular fibrillation, that cause SCA.

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are designed specifically for infrequent use and for use by people whose only training is in CPR and in using AEDs.

Indications for Use

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are indicated for use on patients in cardiac arrest. The patient must be unresponsive (unconscious), not breathing normally, and showing no signs of circulation (for example, no pulse, no coughing, or no movement). With Infant/Child Reduced Energy Defibrillation Electrodes, the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators may be used on children who are up to 8 years old or who weigh up to 25 kg (55 lb).

Contraindications

None known.

Why the Need for Defibrillators

The American Heart Association estimates that, in the USA alone, at least 250,000 people die each year of cardiac arrest. Of these, about 10,000 people might have been saved had they received immediate treatment from a defibrillator.

Sudden cardiac arrest is usually caused by a malfunction in the heart's electrical system. Called ventricular fibrillation, this critical condition prevents the heart from pumping blood throughout the body. Ventricular fibrillation can cause death within seconds.

Defibrillation is a relatively simple procedure that involves placing electrode pads on a victim's exposed chest and delivering an electrical shock to the heart. The externally-delivered shock often restores the heart's electrical system to normal rhythm. Combined with CPR, defibrillation provides the most effective care for victims in cardiac arrest.

Terminology

The following terms appear in this manual.

AED Automated External Defibrillator. A device that evaluates the victim's heart

rhythm and delivers an electrical shock to the heart if a shockable rhythm is

detected.

Cardiac arrest The termination of the heart's pumping action resulting in the lack of a

heartbeat or pulse and breathing.

CPR Cardiopulmonary resuscitation. This involves delivering rescue breathing and

chest compressions to a victim in cardiac arrest.

Defibrillation Delivery of an electrical shock to the heart for the purpose of reversing

ventricular fibrillation.

ECG Electrocardiogram. A composite picture of what is occurring electrically in the

heart.

Fibrillation Chaotic activity of the heart's electrical system. This condition can occur in

the atria or the ventricles. When it occurs in the ventricles, they quiver in a rapid, chaotic manner, preventing them from pumping blood to the body.

Heart attack A nonspecific term referring to the death of heart muscle resulting from

interruption of blood supply, often confused with cardiac arrest.

Impedance Resistance to the flow of electrical current through the body.

Joule The basic unit of energy delivered by a defibrillator.

LED Light emitting diodes.

Myocardial infarction The specific term for what is usually meant by heart attack; death of heart

muscle resulting from an interruption of the blood supply to that area of

myocardium.

Nonshockable rhythm A heart rhythm that is detected by the defibrillator that does not need a

shock, but may need CPR.

Victim In this manual, the person suffering from cardiac arrest.

Responder In this manual, the person giving aid to a victim in cardiac arrest. Used

interchangeably with user.

SAS Physio-Control patented Shock Advisory System™.

Shockable rhythm A heart rhythm that is detected by the defibrillator as requiring a shock, for

example, ventricular fibrillation.

User In this manual, the person giving aid to a victim in cardiac arrest. Used

interchangeably with responder.

Ventricular fibrillation A life-threatening chaotic heart rhythm.

Ventricular tachycardia Rapid heart rhythm originating in the ventricle.

Text Conventions

Throughout this manual, special text characters are used to indicate labels and voice instructions:

Operating control labels: ${\tt CAPITAL}$ LETTERS such as ${\tt ON/OFF}$ and ${\tt SHOCK}.$

Voice instructions: ITALICIZED CAPITAL LETTERS such as STAND CLEAR.

SAFETY INFORMATION

This section provides important information to help you safely operate your defibrillator. Familiarize yourself with all of the terms, warnings, and symbols presented in this section.

Safety Terms

You may encounter the following terms in this manual and while using your defibrillator:

Danger Immediate hazards that will result in serious personal injury or death to the user and/or

the victim.

Warning Hazards or unsafe practices that could result in serious personal injury or death to the

user and/or the victim.

Caution Hazards or unsafe practices that could result in minor personal injury to the user and/or

the victim, product damage, or property damage.

General Warnings and Cautions

WARNINGS!

Shock hazard.

The defibrillator delivers up to 360 joules of electrical energy. Unless properly used as described in this manual, this electrical energy may cause serious injury or death. Do not attempt to operate this defibrillator unless thoroughly familiar with this manual and the function of all controls, indicators, connectors, and accessories.

Do not insert a finger or any object other than the CHARGE-PAK™ into the well of the defibrillator.

Shock or fire hazard.

Do not immerse any portion of this defibrillator in water or other fluids. Avoid spilling any fluids on the defibrillator or its accessories. Do not clean with ketones or other flammable agents. Do not autoclave or sterilize this defibrillator or the accessories unless otherwise specified.

Possible fire or explosion.

Do not store this defibrillator in the presence of flammable gases or in direct contact with flammable material.

Do not use this defibrillator in the presence of flammable gases or anesthetics. Use care when operating this defibrillator close to oxygen sources (such as bag-valve-mask device or ventilator tubing). Turn off gas source or move source away from patient during defibrillation.

Possible device shutdown.

When the attention symbol ♠ first appears on the readiness display, only a small number of shocks and monitoring time are available. Always keep a CHARGE-PAK™ in the defibrillator. Routinely check that the defibrillator is ready for use. Replace the CHARGE-PAK battery charger after each use of the defibrillator.

WARNINGS!

Possible electrical interference with device performance.

Equipment operating in close proximity may emit strong electromagnetic or radio frequency interference (RFI) which could affect the performance of this device. RFI may result in improper device operation, distorted ECG, failure to detect a shockable rhythm, or cessation of pacing. Avoid operating the device near cauterizers, diathermy equipment, cellular phones, or other portable and mobile RF communications equipment. Maintain equipment separation of at least 1.2 m (4 ft) and do not rapidly key EMS radios on and off. Contact a technical support representative if assistance is required.

Possible electrical interference.

Using cables, electrodes, or accessories not specified for use with this device may result in increased emissions or decreased resistance to electromagnetic interference which could affect the performance of this defibrillator or of equipment in close proximity. Use only parts and accessories specified in this manual.

Possible electrical interference.

This device may cause electromagnetic interference (EMI) especially during charge and energy transfers. EMI may affect the performance of equipment operating in close proximity. If possible, verify the effects of defibrillator discharge on other equipment prior to using the defibrillator in an emergency situation.

Possible improper device performance.

Using other manufacturers' cables or electrodes may cause the defibrillator to perform improperly and invalidates the safety agency certification. Use only the parts and accessories specified in this manual.

Using damaged or expired equipment or accessories may cause the defibrillator to perform improperly and may injure the victim or the user.

Safety risk and possible equipment damage.

Monitors, defibrillators, and their accessories (including electrodes and cables) contain ferromagnetic materials. As with all ferromagnetic equipment, these products must not be used in the presence of the high magnetic field created by a Magnetic Resonance Imaging (MRI) device. The high magnetic field created by an MRI device will attract the equipment with a force sufficient to cause death or serious personal injury to persons between the equipment and the MRI device. This magnetic attraction may also damage the equipment. Skin burns will also occur due to heating of electrically conductive materials, such as patient leads and pulse oximeter sensors. Consult the MRI manufacturer for more information.

Possible equipment damage.

This defibrillator may be damaged by mechanical or physical abuse, such as immersion in water or dropping the defibrillator. If the defibrillator has been abused, remove it from use and contact a qualified service technician.

Using damaged or expired equipment or accessories may cause your defibrillator to perform improperly and may injure the victim or the user.

SYMBOLS

The following symbols may appear in this manual and on the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators. For more information about the readiness display symbols, refer to Section 2, Getting Started.

Symbol	Explanation
OK	OK Indicator. The defibrillator is ready for use.
\triangle	Attention. On the readiness display—the internal battery is low. Refer to page 2-4 for more information about the readiness display.
<u>^</u>	Attention. On the CHARGE-PAK battery charger—consult the operating instructions. Refer to page 5-5 for more information about the CHARGE-PAK battery charger.
<u>^</u>	Attention. On the safety warning—consult the operating instructions. Refer to page 1-5 for more information about the warnings and cautions.
<u>^</u>	Attention. On the electrode pads—consult the operating instructions. Refer to page 2-6 for more information about electrode pads.
4	Warning, high voltage.
*	CHARGE-PAK battery charger indicator. On the readiness display—CHARGE-PAK battery charger needs to be replaced.
<i>**</i>	Wrench indicator. There is a condition that prevents or could prevent normal defibrillator operation. Refer to page 5-8 for more information.
	CHARGE-PAK
\uparrow	This end up.
	Fragile/breakable. Handle with care.
**	Protect from water.
①	Power On/Off button.
†	Type BF patient connection.

Symbol

Explanation



Not intended for use on children who are less than eight years of age or who weigh less than 25 kg (55 pounds).



Infant Child Reduced Energy Electrodes are not compatible with QUIK-COMBO defibrillation and therapy cables. To use Infant/Child Electrodes, connect Infant/Child Electrodes directly to the AED.



Not intended for use on adults.



Physio-Control electrodes are latex-free.



Arrow indicates ON/OFF button location.



Symbol denoting a defibrillator and identifies the shock button.



LOT YYWW Lot code.



Do not reuse—single use only.



Use By date shown: yyyy-mm-dd.



Refer to instructions for recycling procedure, page 5-8.



Refer to instructions for disposal procedure, page 5-4.



Do not dispose of this product in the unsorted municipal waste stream. Dispose of this product according to local regulations. See http://recycling.medtronic.com for instructions on the proper disposal of this product.



Mark of conformity according to the European Medical Device Directive 93/42/EEC.



Canadian Standards Association certification for Canada and the United States.

Symbol	Explanation
!USA	For USA audiences only.
760 ————————————————————————————————————	Atmospheric pressure limitation.
595	Relative humidity range 5% to 95%.
50°C - 122°F	Store in a cool, dry location (0° to 50°C or 32° to 122°F).
<1 WEEK - 70°C (158°F) 35°C (95°F) 15°C (69°F) 40°C (-40°F)	Recommended storage temperature: 15° to 35°C (59° to 95°F). Storage at extreme temperatures of -40° or 70°C (-40° or 158°F) is limited to seven days. If storage at these temperatures exceeds one week, the electrode shelf-life will be reduced.
MYYYY	Date of manufacture.
(3)	CAUTION - FIRE HAZARD Do not disassemble, heat above 100°C (212°F), or incinerate battery
\otimes	CAUTION - FIRE HAZARD Do not crush, puncture, or disassemble battery
Rx Only	Federal (U.S.) law restricts this device to sale by, or on the order of, a physician.
MIN	Manufacturer's Item Number.
CAT	Catalog number used for placing orders.
REFor REF	Reorder number.

ABOUT THE LIFEPAK CR PLUS AND LIFEPAK EXPRESS DEFIBRILLATORS

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are designed for indoor and outdoor use. Each has two models available—fully automatic and semiautomatic. After the electrode pads are applied, the fully automatic model evaluates the heart rhythm and, if a shockable rhythm is detected, delivers a shock without any responder assistance. The semiautomatic model evaluates the heart rhythm but requires the responder to press the shock button if a shockable rhythm is detected. Both models have voice instructions that guide the responder through the defibrillation process.

Capabilities and Features

The following paragraphs introduce specific features found in the defibrillators.

Accessories

The defibrillator arrives with an installed CHARGE-PAK battery charger, one preconnected QUIK-PAK™ electrode packet, and the operating instructions. Refer to Section 5 for other accessories.

Automated Operation

Voice instructions guide the responder through the defibrillation process.

The **fully automatic** defibrillator requires no operator interaction beyond placing the electrode pads on the victim. If the defibrillator detects a shockable rhythm, it warns the responder prior to delivering any shock, and then delivers a shock without operator interaction.

The **semiautomatic** defibrillator has a highly visible shock button that the responder must press when the defibrillator issues a shock voice instruction.

Automatic Self-Test

The defibrillator tests itself each week and every time you turn it on. In addition, every month, the defibrillator performs a more extensive self-test. This self-test checks the defibrillator's circuitry to verify that it is ready for use.

Customized Setup

The defibrillator is shipped ready to use with the preprogrammed ADAPTIV biphasic escalating energy protocol. The defibrillator also has several operating settings that can be customized. These include the defibrillator ID, voice prompt volume, energy sequence, and protocol.

Data Management

The defibrillator digitally stores data when it is turned on and the electrode pads are successfully applied to the victim. The stored data includes date and time, ECG data, and the number of shocks. The defibrillator also stores the results of the automatic self-tests.

Stored data can be transferred to a PC by means of a serial infrared link, the IrDA® port. A data transfer and management program running on the PC transfers event and test data from the defibrillator.

Defibrillation Electrodes (Pads)

When applied to the victim, Physio-Control QUIK-PAK™ defibrillation electrodes (pads) work with the defibrillator to monitor the heart rhythm and identify when a shock should be delivered. If victim care is transferred to emergency medical personnel, these electrode pads can be disconnected from the defibrillator and reconnected to other AEDs or defibrillators that are compatible with QUIK-COMBO™ electrodes. For infants or children who are less than eight years of age or weigh less than 25 kg (55 lbs), use Infant/Child Reduced Energy Defibrillation Electrodes. These electrodes reduce the energy delivered by the AED to the victim by one fourth. Keep all electrode pads with the AED.

Defibrillation Waveform

The defibrillation shock, using ADAPTIV Biphasic technology, is delivered in the form of a biphasic truncated exponential (BTE) defibrillation waveform.

Heart Rhythm Analysis

The Physio-Control patented Shock Advisory System™ evaluates the victim's heart rhythm. Refer to Appendix B for further information.

Motion Detection

This patented motion system detects victim or operator motion that could affect the heart rhythm evaluation. Heart rhythm evaluation is interrupted if the defibrillator detects motion.

Power System

Power is provided by a rechargeable internal lithium battery. The internal battery supplies power to operate the defibrillator. A replaceable CHARGE-PAK battery charger provides a trickle charge for the internal battery. It is important to keep a fresh CHARGE-PAK battery charger in the defibrillator, even when the defibrillator is stored.

Readiness Display

This easy-to-read visual display indicates if the defibrillator is ready for use or if it needs attention.

GETTING STARTED

This section provides an orientation to the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators and describes how to prepare the defibrillator for use.

Unpacking and Inspecting Your LIFEPAK CR Plus or LIFEPAK EXPRESS Defibrillator	page 2-2
Where to Locate Your LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator	2-3
Controls, Indicators, and Labels	2-4

UNPACKING AND INSPECTING YOUR LIFEPAK CR PLUS OR LIFEPAK EXPRESS DEFIBRILLATOR

To help ensure the integrity of your defibrillator and to verify that it is complete, perform the initial inspection as follows:

- 1 Remove your defibrillator and examine its outside for signs of damage that may have occurred during shipping.
- 2 Check the remaining contents in the box against the sales order.
- 3 View the OK symbol in the readiness display.

 This indicates your defibrillator is ready for use. If the OK symbol is not visible, contact your local Physio-Control representative.
- 4 Notice the Use By date.
 - The Use By date is located below the readiness display. This date tells you when the electrode packet and battery charger must be replaced.
- 5 Check the defibrillator speaker by performing the following:

Note: This is only a speaker check. **Do not respond to the voice instructions.**

- Press the ON-OFF button to open and turn on your defibrillator.
 Confirm that the voice instructions sound.
- Press and hold the ON-OFF button for approximately 2 seconds to turn off your defibrillator. Three tones will sound.
- 6 Close and latch the lid. Do not reopen the lid unless necessary. Doing so will reduce battery power.

If you have any questions about your defibrillator, please call your local Physio-Control representative.

CAUTION!

After completing an initial inspection, do not open the lid unnecessarily. Each time you open the lid, the defibrillator turns on and internal battery power is reduced. After 30 minutes of cumulative on time, the CHARGE-PAK indicator appears on the readiness display indicating the CHARGE-PAK battery charger and the QUIK-PAK electrode packet should be replaced.

Save the shipping container and inserts in case you need to reship the defibrillator in the future.

WHERE TO LOCATE YOUR LIFEPAK CR PLUS OR LIFEPAK EXPRESS DEFIBRILLATOR

The defibrillator should be easy to reach in a location free of obstacles. This could include a location near existing emergency equipment, such as fire extinguishers and first-aid kits. When considering location, avoid areas that expose the defibrillator to moisture, dust, or extreme temperatures. Recommended storage temperature is 15° to 35° C (59° to 95° F). Storage at higher temperatures will shorten the life of the battery and electrodes.

WARNING!

Possible fire or explosion.

Do not store this defibrillator in the presence of flammable gases or in direct contact with flammable material.


Although the defibrillator and electrodes are designed to withstand environmental temperature fluctuations between -40° C to 70° C (-40° F to 158° F), storage at extreme temperatures of -40° C or 70° C (-40° F or 158° F) is limited to one week. If storage at these temperatures exceeds one week, the electrode shelf-life will be reduced. Refer to Appendix A, page A-3 environmental specifications information. You can place your defibrillator on a stable surface, or you can mount it using the wall mount bracket accessory. Contact your local Physio-Control representative.

CONTROLS, INDICATORS, AND LABELS

This section introduces you to the controls, indicators, and labels on the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators.

Outside Controls, Indicators, and Labels

Controls, indicators, and labels on the outside of the defibrillator are identified in Figure 2-1 and described in Table 2-1.

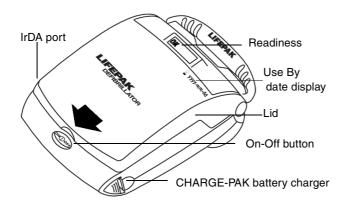

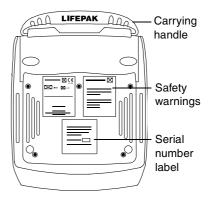


Figure 2-1 Outside Controls, Indicators, and Labels

Table 2-1 Outside Controls, Indicators, and Labels

Description	
There are four indicators that can appear when your defibrillator is turned off. These indicators allow you to determine, just by looking at the defibrillator, whether it's ready for use or needs attention. They include the following:	
OK The OK indicator appears when the defibrillator is turned off and ready for use.	
The CHARGE-PAK indicator appears when the CHARGE-PAK battery charger needs to be replaced or is not installed in the defibrillator. If needed, the defibrillator can be used in an emergency.	
The attention indicator appears when the internal battery is not fully charged. When this indicator first appears, the internal battery can power the defibrillator for a minimum of 6 shocks or 42 minutes.	
The wrench indicator appears when a condition prevents or could prevent the defibrillator from operating normally.	
The top of the defibrillator.	
The ON-OFF button opens the defibrillator lid and turns the defibrillator on. Pressing and holding the button for approximately 2 seconds after the lid is open turns off the defibrillator.	

Table 2-1 Outside Controls, Indicators, and Labels (Continued)

Feature	Description
CHARGE-PAK battery charger	The CHARGE-PAK battery charger delivers a trickle charge to the internal battery. The battery charger can provide a charge for approximately two years, as long as the defibrillator is not used.
IrDA port	The Infrared Data Association defines specifications for infrared wireless communications. The IrDA port provides wireless communications for transferring data from your defibrillator to a PC.
Carrying Handle	The carrying handle is used to transport the defibrillator.
Safety Warnings	Safety warnings provide important information concerning the defibrillator's use and service.
Serial Number Label	The serial number label includes the defibrillator identification number.

Inside Features

The inside features of the defibrillator are designed to make it easy to use during a cardiac arrest event. When you press the ON-OFF button, the lid opens, the defibrillator turns on, and you see the electrode packet and its release handle as shown in Figure 2-2. Table 2-2 describes the inside features.

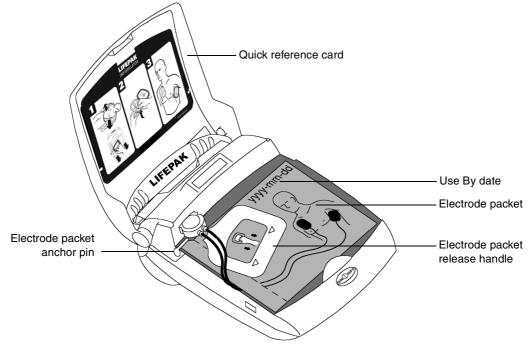


Figure 2-2 Inside Features

Table 2-2 Inside Features

Feature	Description
Quick reference card	This card provides abbreviated graphic directions for using your defibrillator to treat a victim in cardiac arrest.
Use By date	Use By date shown (yyyy-mm-dd) can be viewed through the defibrillator lid when it is closed.
Electrode packet	The QUIK-PAK electrode packet is preconnected to the defibrillator. This packet contains a set of electrode pads.
Electrode packet release handle	When you pull this handle, the electrode packet tears open.
Electrode packet anchor pin	This pin securely positions the electrode packet to the defibrillator.

After you pull the electrode packet release handle and tear open the electrode packet, you will see the features shown in Figure 2-3

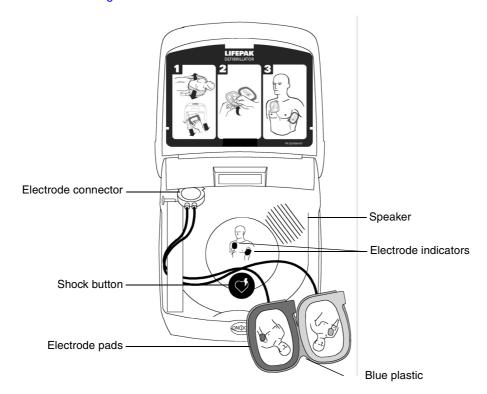


Figure 2-3 Inside Features After Releasing the QUIK-PAK Electrode Packet

Table 2-3 Inside Features After Releasing the QUIK-PAK Electrode Packet

Feature	Description
Speaker	This projects the voice instructions that lead you through the defibrillation process.
Electrode indicators	The electrode indicators flash red until the pads are applied to the victim's exposed chest. When the pads are successfully applied, the indicators turn solid green and the defibrillator can perform an analysis.
	In addition, electrode indicators briefly flash when the defibrillator performs an automatic self-test.
Blue plastic	The plastic liner protects the conductive adhesive gel until the electrode pads are used.
Electrode pads	The electrode pads are applied to the victim's exposed chest; they transfer the defibrillation energy (shock) to the victim. The electrode pads must be removed from the blue plastic before applying them to the victim.
SHOCK button	The SHOCK button is only provided on the semiautomatic model. When pressed, this button delivers a shock to the victim. You cannot deliver a shock to a victim unless the defibrillator instructs you to do so.
Electrode connector	The electrode connector is used to connect the electrode pads to the defibrillator. To aid in victim transport, the connector can be unplugged from the defibrillator and plugged into another AED or defibrillator equipped for QUIK-COMBO electrodes.

USING THE DEFIBRILLATOR

This section provides information and instructions for using the LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator on a victim in cardiac arrest.

Warnings and Cautions	page 3-2
Responding to a Sudden Cardiac Arrest emergency	3-2
Voice Instructions and Tones	3-5
Troubleshooting	3-6

WARNINGS AND CAUTIONS

To help ensure safe use of the defibrillator, completely familiarize yourself with the following warnings and cautions.

WARNINGS!

Incorrect rhythm interpretation.

Performing CPR or otherwise handling or transporting the victim while the defibrillator is evaluating the heart rhythm can cause an incorrect or delayed diagnosis. Keep the victim as still as possible while the defibrillator is attached and do not transport the victim.

Shock hazard.

When instructed "Do not touch patient," "Stand by," or "Everyone clear," remain still, do not touch the defibrillator, patient, electrode pads or any material in contact with the patient. Make sure no one is touching the patient when the defibrillator shocks the patient.

Shock hazard.

To remove an unwanted charge, disconnect the electrode cable from the defibrillator, wait for the defibrillator to automatically remove the charge, or turn off the defibrillator.

Possible fire, burns, and ineffective energy delivery.

During defibrillation, material in contact with the electrode pads can cause electrical sparks, skin burns, and divert important defibrillating energy away from the heart. Place electrode pads so that they adhere to the skin completely. Do not allow the electrode pads to touch each other, medication patches, dressings or any other material on the patient's chest.

During defibrillation, air pockets between the skin and electrode pads can cause skin burns. To help prevent air pockets, make sure electrode pads completely adhere to the skin. Do not use damaged, expired, or dried-out electrode pads.

CAUTION!

Possible equipment damage.

Before using the defibrillator, disconnect all equipment from the patient that is not defibrillator-protected.

RESPONDING TO A SUDDEN CARDIAC ARREST EMERGENCY

If not treated, sudden cardiac arrest (SCA) will cause death. In an SCA situation, it is important to remember to immediately call for help and activate your emergency response system.

Basic Steps for Using the LIFEPAK CR Plus or LIFEPAK EXPRESS Defibrillator

Responding to an SCA emergency using the defibrillator involves these basic steps:




Determine if the victim is in SCA. A person in SCA will not respond when you try to shake him or her.



Check for breathing by listening next to the victim's mouth and looking for chest movement.

Use your defibrillator only if the victim is not responding, not moving, and not breathing normally or not breathing at all. If in doubt, use your defibrillator.

Place your defibrillator near the victim and on the side next to you. Press the ON/OFF button to open the lid and turn on your defibrillator. Remain calm. Your defibrillator will guide you through the defibrillation process.

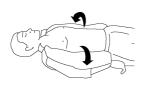

Expose the victim's chest. If the chest is excessively hairy, quickly shave the hair in the area where you will place the pads. If the chest is dirty or wet, wipe the chest clean and dry. If there are medicine patches on the victim's chest, remove them.



Hold down the left side of the electrode packet with one hand and pull the red packet handle down with the other. The electrode packet tears open.

Tear open the packet completely to remove the pads. A small piece of the packet will remain attached to your defibrillator.



Separate the electrode pads, one at a time, from the blue plastic. Use these pads on adults or children 8 years of age or more, who weigh 25 kg (55 pounds) or more. For infants or children who are less than 8 years of age or who weigh less than 25 kg (55 pounds), special electrodes are needed. Refer to page 5-8 for more information.

WARNING!

If you cannot determine a child's age or weight, or if special infant/child electrodes are not available, proceed with the existing electrode pads, and continue on to the next step.



Apply the electrode pads to the victim's bare chest (exactly as shown in the picture on the pads). Be sure to press firmly so that the pads completely adhere to the victim's chest.

Note: Be sure you do not place the electrode pads over an implanted device such as an implanted pacemaker or ICD. An indication of an implant is a protrusion in the chest skin and a scar. If you are in doubt, apply the pads as shown on the labels.



Listen to voice instructions and do not touch the victim unless instructed to do so.



If the defibrillator heart rhythm analysis determines that a shock is needed, the defibrillator will announce PREPARING TO SHOCK, and then instruct you to PRESS FLASHING BUTTON to administer a shock (semiautomatic model) or it will announce PREPARING TO SHOCK, and then automatically administer a shock without requiring further action (fully automatic model).

Do not touch the victim while a shock is delivered.

Regardless of which model you have, continue to follow the voice instructions.



Do not remove the pads or disconnect them from the defibrillator until emergency medical personnel arrive. If the victim starts moving, coughing, or breathing regularly, place the victim in the recovery position (as instructed in CPR training) and keep him or her as still as possible.

What to Do After Emergency Medical Personnel Arrive

When emergency medical personnel arrive, tell them what actions you have taken. Tell them how long the victim has been unconscious, if you delivered shocks, the number of shocks delivered, and if you performed CPR.

Do not worry if you cannot recall precisely what happened. Your defibrillator makes a digital recording of heart rhythms and shocks that can be transferred to a computer at a later time. Refer to Section 4 for information on transferring victim data.

Without removing the electrode pads from the victim, emergency medical personnel can disconnect the electrode pads from the defibrillator and reconnect them to another defibrillator or AED that has a compatible QUIK-COMBO cable.

To disconnect the electrode pads:

- 1 Pull the electrode cable straight out from the defibrillator.
- 2 Remove the electrode packet anchor pin from the slot in the defibrillator.
- 3 Press the ON-OFF button and close the lid to turn off the defibrillator.

What to Do After Using Your Defibrillator

After you use your defibrillator to respond to an SCA emergency, complete the following tasks:

- 1 If the defibrillator is turned on, press and hold the ON-OFF button for approximately 2 seconds to turn it off
- 2 Clean the defibrillator and its accessories according to the instructions provided in Table 5-1, page 5-3. Use only the cleaning agents identified in Table 5-1.
- 3 Transfer data, if desired.
- 4 Replace the CHARGE-PAK battery charger. (Refer to page 5-5.)
- 5 Install a new QUIK-PAK electrode packet. (Refer to page 5-7.)
- 6 Close the lid and verify that the OK symbol appears in the readiness display, indicating that the defibrillator is ready for use. If the attention symbol appears after you replace the battery charger, the internal battery needs additional time to reach an adequate charge capacity.
- 7 Dispose of the used electrode pads, any unused spare electrode pads, and the battery charger. (Refer to Recycling Information, page 5-8.)

VOICE INSTRUCTIONS AND TONES

Defibrillator voice instructions provide clear, step-by-step instructions for responding to a victim in cardiac arrest. In addition, your defibrillator may emit sounds that alert you to the actions that the defibrillator is performing.

Note: A few seconds may pass between voice instructions and tones. Always wait for further instructions before taking action.

Note: Some voice instructions will repeat during the defibrillation process.

TROUBLESHOOTING

This section explains problems you may encounter while using the LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator. For information about keeping your defibrillator in a state of readiness, refer to Section 5.

Table 3-1 Troubleshooting During Victim Use

Problem	Possible Cause	What To Do
CHECK PADS FOR GOOD CONTACT or CHECK CONNECTOR voice	Inadequate connection to the defibrillator	Be sure that the connector is completely inserted.
instructions are heard	Electrode pads are not properly adhered to the victim	 Firmly press the pads on the victim's skin. Clean, shave, and dry the victim's skin before placing pads on skin.
	Electrode pads are dry, damaged, or have passed the expiration date	Replace the pads.
	Electrode pads are not removed from the blue plastic	 Remove the pads from the blue plastic and apply them to the victim's chest.
Defibrillator cannot deliver the required shock	Defibrillator internal battery power is low	Administer CPR if the victim is not responding, not breathing normally, or not moving.
Voice instructions sound faint or distorted	Defibrillator internal battery power is low	Administer CPR if the victim is not responding, not breathing normally, or not moving.
MOTION DETECTED and STOP MOTION voice	Victim movement because of location	Move the victim to a stable location, if possible.
instructions are heard	Victim movement because of breathing	Stop CPR during analysis.Check victim for normal breathing.
	Vehicle motion	 Stop vehicle during analysis, if possible.
	Electrical/radio frequency interference	 Move communication or other suspected devices away from the defibrillator when possible.

Table 3-1 Troubleshooting During Victim Use (Continued)

Problem	Possible Cause	What To Do
Defibrillator does not deliver voice instructions or beeping tones after you open it (turn it on)	Depleted internal battery	 Administer CPR if the victim is not responding, not breathing normally, or not moving.
(tam it on)		 Replace the CHARGE-PAK battery charger as soon as possible. After the OK symbol appears on the readiness display, return the defibrillator to service.
		 Contact authorized service personnel.
	Speaker system failure	 Administer CPR if the victim is not responding, not breathing normally, or not moving.
		 Contact authorized service personnel.
The readiness display is blank	The defibrillator has been turned on	Normal condition when the defibrillator is in use.
	Operating temperature is too low or too high	 Operate the defibrillator within 0° to 50°C (32° to 122°F).
	LCD not operating properly	 Contact authorized service personnel.

DATA STORAGE

This section describes the data that the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators store when used in an SCA event.

This manual does not provide instructions for how to hand off defibrillation data to the emergency medical system or hospital personnel. Because this process varies from area to area, check with the emergency medical system administration for information and directions.

Overview of Data Storage

page 4-2

OVERVIEW OF DATA STORAGE

Each time you use the LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator, it digitally saves data about the victim that can be transferred to a PC. This data can be provided to emergency medical personnel or hospital personnel to aid in case review for quality control, training, and research purposes. We recommend that you become familiar with their local requirements for reporting a use of the defibrillator and for providing use data. For assistance in retrieving data from the defibrillator, contact your local Physio-Control sales or service representative.

Data Stored by Your Defibrillator

Whenever you turn on your defibrillator and connect it to a victim, it automatically stores data about the victim. When this data is transferred to a data management system for review, three reports are available: an Event Log, Continuous ECG, and a Summary. Table 4-1 describes these reports.

Table 4-1 Patient Reports

Report Type	Description
Event Log	A chronological log of all events. An event is a condition noted by the defibrillator. Events are listed on page 4-3.
Continuous ECG	Twenty minutes of the victim's ECG rhythm beginning when the victim is connected to the defibrillator and ending when the defibrillator is turned off.
Summary	Combines the Event Log and a sampling of continuous ECG rhythms associated with certain events.

Your defibrillator can store up to two records: one for the current victim and one for the previous victim. When you use your defibrillator, it is important that you transfer this data as soon as possible after use to free up storage.

The Complete Record for the current victim includes the Continuous ECG and Event Log. If you treat a second victim, the first victim's Complete Record will be reformatted into a Summary report. If you treat a third victim, all of the first victim's data will be deleted and the second victim's Complete Record will be reformatted into a Summary Report. Refer to Table 4-2.

Table 4-2 Defibrillator Patient Records

	Complete Record	Summary
Current Victim	~	V
Previous Victim	Ø	•

If you turn your defibrillator on and off without attaching electrode pads to a victim, the defibrillator will not create a new record and the records in the defibrillator will not be altered. The defibrillator deletes previous data only after it is connected to a new victim.

After you transfer data records to a PC, the defibrillator will disallow repeat transmissions. However, service personnel may access device records, if necessary.

Test and Service Data

Your defibrillator stores a test log consisting of the most recent auto-tests, power cycles, and CHARGE-PAK battery charger replacements. The test log lists the test results and any errors detected. The test log data is available only to service personnel or to users through the data management system.

Event and Test Log

Table 4-3 lists the types of events that may be annotated on event and test log reports.

Table 4-3 Event and Test Log Reports

Event Log	Test Log
Power On	Self-Test Power On
Connect Electrodes	Self-Test Pass/Fail
Patient Connected	User Power On
Initial Rhythm*	CHARGE-PAK Changed
Analysis X*	CHARGE-PAK
Shock Advised	Fault Log
Charge Complete	
SHOCK X-XXXJ*	
Shock X Abnormal	
No Shock Advised	
CPR Prompt	
Stop CPR Prompt	
Check Patient*	
Charge Removed	
Low Battery	
Motion	
Analysis Stopped*	
Out of Event Memory	
Out of Waveform Memory	
Power Off	

^{*}These events include ECG samples in the Summary Report.

CARING FOR YOUR DEFIBRILLATOR

This section explains how to help keep your LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator in good working condition. Cared for properly, your defibrillator is built to give you many years of service.

Maintaining a State of Readiness	page 5-2
Storing Your Defibrillator	5-2
Cleaning Your Defibrillator	5-2
Replacing the CHARGE-PAK Battery Charger and the QUIK-PAK Electrode Packet	5-4
Obtaining Authorized Service	5-8
Recycling Information	5-8
Supplies, Accessories, and Training Tools	5-8
Warranty Information	5-8

MAINTAINING A STATE OF READINESS

Your LIFEPAK CR Plus or LIFEPAK EXPRESS defibrillator does not require routine maintenance. It performs an automatic self-test once a week and every time you turn it on. The electrode indicators briefly flash during the test. If the automatic self-test detects a condition that requires attention, the OK symbol in the readiness display will fade and either the CHARGE-PAK symbol, the ATTENTION symbol, or the WRENCH symbol will appear, depending on the type of condition detected.

On a regular basis, you should do the following:

- Check to make sure that the OK symbol is visible in the readiness display.
- Check the Use By date on the electrode packet (visible through the defibrillator lid in the upper right corner) and all other electrode packets. If the date has passed, replace the electrode packet and the battery charger.
- Check other emergency supplies that may be stored with your defibrillator.

When establishing your local inspection schedule, consider how often your defibrillator will be used and how familiar the operators are with using a defibrillator. For example, if the defibrillator is used only rarely, monthly inspections may be appropriate. An inspection checklist is provided in Appendix C.

STORING YOUR DEFIBRILLATOR

Always store your defibrillator according to the temperature ranges recommended in Appendix A, Specifications.

WARNING!

Possible fire or explosion.

Do not store this defibrillator in the presence of flammable gases or in direct contact with flammable material.

CLEANING YOUR DEFIBRILLATOR

CAUTION!

Possible equipment damage.

Do not clean any part of the defibrillator or accessories with bleach, bleach dilution, or phenolic compounds. Do not use abrasive or flammable cleaning agents. Do not steam, autoclave, or gas-sterilize the defibrillator or its accessories.

Table 5-1 Cleaning Methods

Item	Cleaning Method	Cleaning Agent
Exterior case, readiness	Clean with damp sponge or cloth	Nonabrasive soap and water
display, and crevices		Quaternary ammonium compounds
		Rubbing (isopropyl) alcohol
		Peroxide (peracetic acid) solutions
CHARGE-PAK battery charger	None	None, dispose of/recycle after use
Electrode Pads	None, do not remove electrode pads from the packet	None, dispose of/recycle after use
Carrying case	Wipe with damp cloth or sponge	Water
Quick Reference Card	Wipe with damp cloth or sponge	Water

REPLACING THE CHARGE-PAK BATTERY CHARGER AND THE QUIK-PAK ELECTRODE PACKET

The CHARGE-PAK battery charger is a replaceable, nonrechargeable battery cell that charges your defibrillator's internal battery. The internal battery provides the power for your defibrillator. To prevent damage to the internal battery, always keep the battery charger in place, including during storage or shipping.

The QUIK-PAK electrode packet contains the pads that transfer the defibrillation energy to the victim. The packet should remain connected to the defibrillator and unopened until required for an SCA emergency. QUIK-PAK electrode pads are not reusable.

When installed, these two accessories allow your defibrillator to stand by for use for approximately two years. The electrode packet Use By date is programmed into the battery charger. When the date is reached, the CHARGE-PAK symbol appears in the readiness display, indicating both the battery charger and electrode packet need to be replaced.

Note: The CHARGE-PAK battery charger and QUIK-PAK electrodes are a set and have the same expiration date. Always replace the CHARGE-PAK battery charger and QUIK-PAK electrodes at the same time to keep the replacement cycles in sync.

Use the Physio-Control replacement kit to replace the CHARGE-PAK battery charger and the QUIK-PAK electrode packet as follows:

- · After using the defibrillator
- If the CHARGE-PAK symbol appears in the readiness display
- · When the Use By date is reached or passed

The replacement kit includes a CHARGE-PAK battery charger, one or two QUIK-PAK electrode packets, and a CHARGE-PAK battery discharger. The discharger depletes a used battery charger so that it is ready for recycling or disposal.

Note: Always keep the CHARGE-PAK battery charger with its original defibrillator. When a CHARGE-PAK battery charger is inserted into a defibrillator, it becomes internally linked to that device. Do not move the CHARGE-PAK to another defibrillator. If the CHARGE-PAK is removed from its original defibrillator and inserted into a second device, it will not charge and the CHARGE-PAK indicator will illuminate on the readiness display. If this occurs, discard the existing CHARGE-PAK battery charger and QUIK-PAK electrodes. Obtain a new replacement kit and replace both CHARGE-PAK battery charger and QUIK-PAK electrodes at the same time.

Follow instructions provided in the replacement kit for battery charger and electrode packet recycling/disposal. Contact your local Physio-Control representative.


WARNING!

Possible explosion or fire.

The CHARGE-PAK battery charger is not rechargeable. Do not attempt to recharge, open, crush, or burn the battery, or it may explode or catch fire.

Replacing the CHARGE-PAK Battery Charger

To replace the CHARGE-PAK Battery Charger:

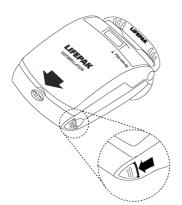

1 Press the release latch (in the direction of the arrow) to remove the used battery charger. The battery charger springs outward from the defibrillator.



- 2 Insert the new battery charger into your defibrillator and push until you hear it click into position.
- 3 Confirm that the symbol disappears and that the **OK** symbol appears in the readiness display.

Note: If the symbol appears after you replace the battery charger, the internal battery is very low and needs time to charge. It may take up to three days if you had the defibrillator on for a long time or if you delivered many shocks. The **OK** symbol appears when the internal battery is charged.

Remember: If the defibrillator is needed for an emergency, attempt to use it even if the symbol is visible.

CAUTION!

Keep the defibrillator at temperatures between $0^{\circ} - 50^{\circ}$ C ($32^{\circ} - 122^{\circ}$ F) while the new battery charger charges the internal battery. The internal battery may not charge effectively at lower temperatures. Temperatures exceeding 50° C (122° F) for longer than seven days can permanently damage the internal battery.

Caring for Your Defibrillator

To discharge and dispose of a used battery charger:



1 Insert the discharger into the used battery charger.



Note: Do not attempt to remove the discharger once it is in place.



2 Let the discharger fully deplete the battery charger. Wait at least 9 days.



3 Place the used battery charger in the trash or recycle it.

Replacing the QUIK-PAK Electrode Packet


To replace the QUIK-PAK Electrode Packet:









2 Press and hold down the ON/OFF button for 2 seconds to turn off the defibrillator and save battery power.

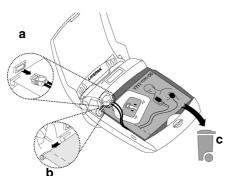

- 3 Remove the outdated or used electrode packet:
 - a Unplug the electrode connector from the connector receptacle.
 - b Slide the anchor pin from the slot.
 - Discard the outdated or used electrode packet according to local regulations.

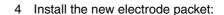

- a Slide the anchor pin into the slot.
- b Plug the electrode connector into the receptacle.
- c Ensure that the new electrode packet is centered on the defibrillator and is tucked behind the lip before closing the lid.
- d Close the lid. Confirm that the packet Use By date is visible through the upper right-hand corner of the lid.

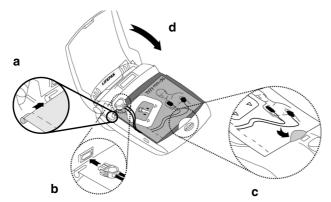

OBTAINING AUTHORIZED SERVICE

WARNING!

Shock hazard.

Do not disassemble the defibrillator. It contains no operator-serviceable components and dangerous high voltages may be present. Contact authorized service personnel for repair.

The WRENCH indicator appears on the readiness display if the defibrillator requires service. Contact authorized service personnel only. In the USA, contact Physio-Control at 1.800.442.1142. Outside the USA, contact your local Physio-Control representative. Be prepared to provide the following information:

- Model number and MIN number (part number)
- Serial number
- · Description of the problem based on your observations

RECYCLING INFORMATION

Recycle the defibrillator and its accessories at the end of their useful lives.

Recycling Assistance

Items should be recycled according to national and local regulations. Contact your local Physio-Control representative for assistance or refer to http://recycling.medtronic.com for instructions on disposing of this product.

Preparation

Items should be clean and contaminant-free prior to being recycled.

Recycling of Disposable Electrodes

After using disposable electrodes, follow your local clinical procedures for recycling.

Packaging

Packaging should be recycled according to national and local regulations.

SUPPLIES, ACCESSORIES, AND TRAINING TOOLS

If possible, we recommend that you have an extra QUIK-PAK Electrode Packet and CHARGE-PAK Battery Charger on hand. In addition, there are other useful accessories available. For example, you can buy Infant/Child Reduced Energy Defibrillation Electrodes for use with your defibrillator on children who are less than 8 years of age or weigh less than 25 kg (55 lb). Contact your local Physio-Control representative.

WARRANTY INFORMATION

Refer to the warranty statement included in the accessory kit shipped with this product. For duplicate copies, contact your local Physio-Control representative. In the USA, call 1.800.442.1142. Outside the USA, contact your local Physio-Control representative.

DEFIBRILLATOR OPERATING SETTINGS

This section introduces the operating settings that are adjustable on the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators.

Operating Settings and Setup Configuration

page 6-2

OPERATING SETTINGS AND SETUP CONFIGURATION

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators have 12 operating settings that affect how it functions. The settings range from defining the time and date used by the defibrillator, to the energy sequence and protocol of the shocks delivered. This group of operating settings are the defibrillator's setup configuration.

Each operating setting is preset to a default value based on clinical guidelines and does not need to be changed unless your medical director desires performance characteristics that are different from the default setting.

Table 6-1 identifies operating settings in the setup configuration, describes each setting, including possible options, and identifies the preset defaults.

Table 6-1 Operating Settings

Operating Settings		Default Setting	
Operating Settings	Description	Default Setting	
Device ID	Device ID is a unique identifier (ID) assigned to each defibrillator used for tracking the location of defibrillators. If you have multiple defibrillators at your location, you should give each a unique name, such as Building 1, Building 2, and so on. When you transfer event data from the defibrillator to a PC, the defibrillator ID will be included in the transferred data.	Serial Number	
Energy Sequence	Energy sequence defines the energy levels used	Level 1 – 200 joules	
	by the defibrillator.	Level 2 – 300 joules	
	The energy level choices are: 150, 175, 200, 225, 250, 275, 300, 325, 360.	Level 3 – 360 joules	
Energy Protocol	Energy protocol determines how the defibrillator delivers successive shocks. There are two options for this setting: flexible or fixed.	Flexible s	
	Flexible sequence means the energy delivered for a shock increases only if an analysis immediately following a shock results in another SHOCK ADVISED decision. For example, if the defibrillator energy sequence is set up as 200, 300, 360, flexible sequence means that the energy delivered for the first shock is 200 joules. If the arrhythmia is terminated by shock 1 and the next analysis results in a NO SHOCK ADVISED decision, the energy will not increase for the next shock. However, if the arrhythmia is not terminated by shock 1 and the next analysis results in a SHOCK ADVISED decision, the energy will increase to 300 joules, and so on. Fixed sequence means that the energy delivered after the first shock of 200 joules increases from 200 to 300, and then to 360 joules, regardless of the post-shock ECG rhythm and subsequent analysis.		

Defibrillator Operating Settings

Table 6-1 Operating Settings (Continued)

Operating Settings	Description	Default Setting
CPR Time 1 CPR Time 2	The CPR Time 1 and CPR Time 2 settings define the time interval for performing CPR after a shock or after a no shock advised decision. The choices for CPR Time 1 and CPR Time 2 are	CPR Time 1–120 sec CPR Time 2–120 sec
	15, 30, 45, 60, 90, 120, and 180 seconds.	
Device Date Device Time	Device date and device time sets the current date and time.	Pacific Standard Time and Date
Turn-On Prompt	This setting determines whether you hear only tones, or tones and the voice prompt, CALL FOR HELP NOW, when you open the lid and the defibrillator turns on. Choices are Voice and Tones.	Voice
Voice Prompt Volume	This setting sets the voice prompt volume to MEDIUM or HIGH.	High
Pulse Prompt	The pulse prompt setting determines whether the defibrillator prompts you to check the victim's pulse (appropriate for medically trained users) or to check the victim for signs of circulation such as breathing and movement (appropriate for lay users). Choices are Check Pulse, Check Breathing, or Check Circulation.	Per customer order
Stack Shocks	When set to OFF, the Stack Shocks option eliminates the analysis after each shock and inserts prompting for CPR after each (a single) shock. This eliminates the three-shock stack. CPR is prompted regardless of the ECG rhythm after the shock. The CPR time following the shock is determined by the CPR Time 1 setting selected. Choices for Stack Shocks option are ON or OFF. When set to ON, an analysis will occur after shocks, and up to three shocks in a row may be delivered (three-shock stack).	Off
Pulse Check	When set to Never, the Pulse Check option removes all prompting for pulse checks. The other Pulse Check settings available allow pulse checks only after every No Shock Advised (NSA) decision, after the second NSA decision and thereafter, or Always (after shocks, NSA and CPR).	Never

Defibrillator Operating Settings

Table 6-1 Operating Settings (Continued)

Operating Settings	Description	Default Setting
Motion Detection	The motion detection setting is used to determine if motion detection is active or not active during analysis.	On
	When motion detection is On, the defibrillator stops analysis for up to 10 seconds if it detects any victim motion. The defibrillator notifies the responder of the problem. The defibrillator will resume analysis after 10 seconds, even if motion is still present. When motion detection is Off, analysis is not inhibited, regardless of any victim motion.	
Time Zone	This setting sets the time zone for the location of the defibrillator. Choices are 74 time zones with universal time code (UTC).	None

APPENDIX A SPECIFICATIONS

SPECIFICATIONS

All specifications are at 20°C (68°F) unless otherwise stated.

Defibrillator

Waveform: Biphasic Truncated Exponential, with voltage and duration compensation for

victim impedance. See Figure A-1.

Output Energy Sequence:

Multiple levels, configurable from 150 joules to 360 joules.

The following specifications apply from 25 to 200 ohms. Voltage compensation is limited to the voltage that would result in delivery of 360 joules into 50 ohms.

Output Energy $\pm 10\%$ into 50 ohms Accuracy: $\pm 15\%$ into 25 to 100 ohms

Waveform Parameters:

Impedance (Ω)	Phase 1 Dui	ration (ms)	Phase 2 Du	ration (ms)	Tilt (%)
	Min	Max	Min	Max	Nominal
25	5.1	6.0	3.4	4.0	77.5
50	6.8	7.9	4.5	5.3	65.8
100	8.7	10.6	5.8	7.1	52.6
125	9.5	11.2	6.3	7.4	47.8

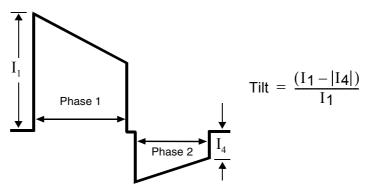


Figure A-1 Biphasic Waveform

Shock Advisory System: An ECG analysis system that advises whether a shock is appropriate, meets

rhythm recognition criteria specified in the American Association of Medical

Instrumentation standard DF39.

The device allows a defibrillation shock only if the Shock Advisory System

advises defibrillation.

Device Capacity: LIFEPAK CR Plus defibrillator — Thirty (30) full discharges or 210 minutes

of "ON time" with a fully charged device.

LIFEPAK EXPRESS defibrillator — Twenty (20) full discharges or

140 minutes of "ON time" with a fully charged device.

Shock Charge Time: Charge times with a fully charged device:

200 joules in less than 9 seconds360 joules in less than 15 seconds

Charge time after 15 discharges delivered from a fully charged device:

360 joules in less than 15 seconds

System Recharge

Times:

Recharge times with a fully discharged device:

Able to deliver 6 shocks or provide 42 minutes of operating time after 24 hours of recharge and 20 shocks or 140 minutes of operating time after 72 hours of recharge time with a new CHARGE-PAK at temperatures above

15°C (59°F).

Controls: LID RELEASE/ON-OFF button—Controls device power.

SHOCK button (semiautomatic version)—Delivers defibrillation energy.

After electrodes are attached to a victim, an automatic version of the device

delivers a shock, if appropriate, not requiring operator intervention.

Electrical Protection: Input protected against high voltage defibrillator pulses per

IEC60601-1/EN60601-1. See Figure A-2.

┤★

Figure A-2 Defibrillation-protected, Type BF Patient Connection

Safety Classification: Internally powered equipment. IEC60601-1/EN60601-1.

User Interface

User Interface: The user interface includes voice instructions, audible tones, and graphical

prompts.

Readiness Display The readiness display shows the device status.

OK







OK Indicator: Indicates OK when the last self test was completed successfully. When the

OK indicator is visible, all other indicators are not visible.

The OK indicator is not displayed during device operation.

CHARGE-PAK

Indicator:

When displayed, replace the CHARGE-PAK battery charger.

Attention Indicator: When first displayed, at least 6 discharges or 42 minutes of operating time

will remain.

Service Indicator: Service required when displayed.

Environmental

Note: All performance specifications defined assume that the unit has been stored (two hours minimum) at the operating temperature prior to operation.

Operating Temperature: 0 to 50° C (32 to 122° F)

Storage Temperature: -40 to 70° C (-40 to 158° F) with CHARGE-PAK and electrodes, maximum

exposure time limited to one week.

Atmospheric Pressure: 760 mmHg to 429 mmHg, 0 to 15,000 feet above sea level.

Relative Humidity: 5 to 95% (non-condensing)

Water Resistance: IEC60529/EN60529 IPX4 "Splash proof" with electrodes connected and

CHARGE-PAK installed.

Shock: MIL-STD-810E, Method 516.4, Procedure 1, (40g, 6 to 9 msec pulse, ½ sine

each axis).

Vibration: MIL-STD-810E, Method 514.4, Helicopter – category 6 (3.75 g rms) and

Ground Mobile - category 8 (2.85 g rms).

EMC: Refer to Appendix D for EMC information as defined in IEC 60601-1-2.

Physical Characteristics

Height: 10.7 cm (4.2 in)
Width: 20.3 cm (8.0 in)

Depth: 24.1 cm (9.5 in), excluding handle

Weight: 2.0 kg (4.5 lb) with CHARGE-PAK and electrodes

Accessories

CHARGE-PAK Battery Charger

Type: Li/SO₂Cl₂ Lithium Sulfuryl Chloride, 11.7V, 1.4 amp-hours.

Replacement: Replace after each use or when CHARGE-PAK indicator is visible (typically

after 2 years)

Weight: 80.5 grams (0.18 lb)

QUIK-PAK Electrode Pads

Pads: Pacing/defibrillation/ECG electrodes.

Pads Packaging: User-intuitive, rapid-release QUIK-PAK electrodes allow the electrode pads

to be pre-connected to the device and protected under a top cover.

Pads Shelf Life: Two years typical. Electrode Shape: Oval-rectangular.

Electrode Size: 11.2 cm (4.4 in) × 18.5 cm (7.3 in)

Lead Wire: 1.067 meters (3.5 feet)

Conductive Adhesive

Gel Contact Area: 82 cm² (12.8 in²)

Maximum Adhesion

Time: 24 hours

Maximum ECG

Monitoring Time: 24 hours

Maximum Number of

Defibrillation Pulses: 50 at 360 joules

Maximum Pacing

Duration: Up to 12 hours

ECG is received from disposable defibrillation electrodes, standard placement (anterior-lateral), or anterior/posterior placement.

Data Storage

Memory Type: Internal digital memory.

ECG Storage: Dual victim Data Storage.

Minimum 20 minutes of ECG stored for the current victim.

Summarized data stored for the previous victim.

Report Types: Continuous ECG—A continuous ECG report for the victim.

Summary—A summary of critical resuscitation events and ECG waveform

segments associated with these events.

Event Log report—A report of time stamped markers, which reflect operator

and device activity.

Test Log report—A device self test activity report.

Capacity: Minimum 200 time stamped Event Log markers.

Communications: Wireless transfer to a personal computer.

Data Review: Physio-Control provides an array of tools to meet customer needs for data

viewing and analysis.

CLINICAL SUMMARY: DEFIBRILLATION OF VENTRICULAR FIBRILLATION AND VENTRICULAR TACHYCARDIA

Background

Physio-Control conducted a multi-centered, prospective, randomized and blinded clinical trial of biphasic truncated exponential (BTE) shocks and conventional monophasic damped sine wave (MDS) shocks. Specifically, the equivalence of 200 J and 130 J BTE shocks to 200 J MDS shocks¹ was tested.

Methods

Ventricular fibrillation (VF) was induced in 115 patients during evaluation of implantable cardioverter defibrillator function and 39 patients during electrophysiologic evaluation of ventricular arrhythmias. After 19±10 seconds of VF, a customized defibrillator delivered an automatically randomized shock. Efficacy was based on success of this shock. To demonstrate equivalence of test shocks to control shocks, the 95% upper confidence limit of the difference in efficacy (95UCLD), control minus test, was required to be less than 10%.

Results

Ventricular Fibrillation

The efficacy of the 200 J BTE shocks was demonstrated to be at least equivalent to the efficacy of 200 J MDS shocks (95UCLD=2%). The difference is success rates of 200 J MDS minus 200 J BTE shocks was -10% (exact 95% confidence interval from -27% to 4%). The 130 J BTE shocks were not demonstrated equivalent to 200 J MDS shocks (95UCLD=22%). However, neither was their efficacy significantly lower than that of the 200 J MDS shocks (statistical power limited by small sample sizes). For all shock types, hemodynamic parameters (oxygen saturation and systolic and diastolic blood pressure) were at or near their pre-induction levels by 30 seconds after successful shocks.

Shock	Ventricular Fibrillation 1st Shock Success	Exact 95% Confidence Interval
200 J MDS	61/68 (90%)	80 to 96%
200 J BTE	39/39 (100%)	91 to 100%
130 J BTE	39/47 (83%)	69 to 92%

Ventricular Tachycardia

Seventy-two episodes of ventricular tachycardia (VT), induced in 62 patients, were treated with randomized shocks. High rates of conversion were observed with biphasic and monophasic shocks. Sample sizes were too small to statistically determine the relationship between success rates of the waveforms tested.

¹S.L. Higgins et al., "A comparison of biphasic and monophasic shocks for external defibrillation," *Prehospital Emergency Care*, 2000, 4(4):305-13.

Shock	Ventricular Fibrillation 1 st Shock Success	Exact 95% Confidence Interval
200 J MDS	26/28 (93%)	77 to 99%
200 J BTE	22/23 (96%)	78 to 100%
130 J MDS	20/21 (95%)	77 to 100%

Conclusions

In this double-blinded study, the efficacy of the 200 J BTE shocks was demonstrated to be at least equivalent to the efficacy of 200 J MDS shocks for defibrillation of short duration, electrically-induced VF. However, the comparison of efficacy of 130 J biphasic and 200 J monophasic shocks for VF was inconclusive. All waveforms tested provided a high rate of termination of VT. The VT sample sizes were too small to statistically determine the relationship between VT success rates of the waveforms tested.

Compared to conventional shocks for VF, we found no positive or negative effect of biphasic shocks for VF on hemodynamic parameters following the defibrillating shock. It is possible that, compared to 200 J monophasic shocks, 200 J biphasic shocks will in some cases enable earlier termination of VF. Therefore, we conclude that biphasic shocks for VF delivered at conventional energy levels have the potential to improve outcome in resuscitation of victims with sudden cardiac arrest.

APPENDIX B SHOCK ADVISORY SYSTEM

OVERVIEW OF THE SHOCK ADVISORY SYSTEM

The Shock Advisory System (SAS) is an ECG analysis system built into the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators that advises the operator if it detects a shockable or nonshockable rhythm. This system makes it possible for individuals not trained to interpret ECG rhythms to provide potentially-lifesaving therapy to victims of ventricular fibrillation or pulseless ventricular tachycardia. The SAS contains the following features:

- Electrode contact determination
- · Automated interpretation of the ECG
- · Operator control of shock therapy
- · Motion detection

Electrode Contact Determination

The victim's transthoracic impedance is measured through the defibrillation electrodes. If the baseline impedance is higher than a maximum limit, it is determined that the electrodes are not in sufficient contact with the victim or not properly connected to the defibrillator. ECG analysis and shock delivery are inhibited. The operator is advised to connect electrodes any time electrode contact is inadequate.

Automated Interpretation of the ECG

The Shock Advisory System is designed to recommend a shock if it detects the following:

- Ventricular fibrillation with a peak-to-peak amplitude of at least 0.08 mV
- Ventricular tachycardia defined as having a heart rate of at least 120 beats per minute, QRS width of at least 0.16 seconds, and no apparent P waves.

Pacemaker pulses may prevent advisement of an appropriate shock, regardless of the victim's underlying rhythm. The SAS is designed to recommend no shock for all other ECG rhythms including pulseless electrical activity, idioventricular rhythms, bradycardia, supraventricular tachycardias, and normal sinus rhythms.

ECG analysis is performed on consecutive 2.7 second segments of ECG. The analysis of two out of three segments must agree before a decision (SHOCK ADVISED or NO SHOCK ADVISED) is made.

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillator SAS performance for adult, pacemaker and pediatric ECGs is summarized in the following table.

Table B-1 LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillator SAS Performance Table for Adult ECGs

Rhythm Class	ECG Test ¹ Sample Size	Performance Goal ^{2, 3}	Observed Performance Sensitivity or Specificity [LCL] ⁴
Shockable: coarse VF	168	>90% sensitivity	100.0% [98.6%]
Shockable: shockable VT	65	>75% sensitivity	84.6% [77.3%]
Nonshockable: NSR	144	>99% specificity for NSR (AHA)	100.0% [98.4%]
Nonshockable: asystole	43	>95% specificity	100.0% [94.8%]

Table B-1 LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillator SAS Performance Table for Adult ECGs (Continued)

Rhythm Class	ECG Test ¹ Sample Size	Performance Goal ^{2, 3}	Observed Performance Sensitivity or Specificity [LCL] ⁴
Nonshockable: all other rhythms	531	>95% specificity	95.9% [94.5%]
Intermediate: fine VF	29	Report only	96.6% [87.2%] sensitivity

¹ From Physio-Control ECG database. Each sample is run 10 times asynchronously.

VF = ventricular fibrillation

VT= ventricular tachycardia

NSR = normal sinus rhythm

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators were also tested using paced rhythms recorded at high-fidelity from victims with implanted pacemakers. The high-fidelity pacemaker spikes were also added to samples of ventricular fibrillation to test the defibrillator's ability to reach a shock decision in the case of ventricular fibrillation with an implanted, active pacemaker. The results are summarized in the following table.

Table B-2 LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillator SAS Performance with Active Pacemakers

Rhythm Class	ECG Test Sample Size	Performance Goal	Observed Performance
Shockable: Coarse VF	35	>90% sensitivity	91.4% [81.9%]
Nonshockable: Paced rhythms	35	>95% specificity	100.0% [93.6%]

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators were also tested using ECGs acquired from hospitalized pediatric victims ranging in age from < 1 day old to 17 years old. The results are summarized in the following table.

 Table B-3
 LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillator SAS Performance Table for Pediatric ECGs

Rhythm Class	ECG Test ¹ Sample Size	Performance Goal ²	Observed Performance Sensitivity or Specificity [LCL] ³
Shockable: coarse VF	90	>90% sensitivity	100.0% [97.5%]
Shockable: shockable VT	11	>75% sensitivity	54.5% [31.8%]
Nonshockable: NSR	424	>99% specificity	100.0% [99.5%]
Nonshockable: asystole	95	>95% specificity	100.0% [97.6%]

² Association for the Advancement of Medical Instrumentation. DF39-1993 Standard for Automatic External Defibrillators and Remote-Control Defibrillators. Arlington, VA: AAMI;1993.

³ Automatic External Defibrillators for Public Access Defibrillation: Recommendations for Specifying and Reporting Arrhythmia Analysis Algorithm Performance, Incorporating New Waveforms, and Enhancing Safety. AHA Task Force on Automatic External Defibrillation, Subcommittee on AED Safety and Efficacy. Circulation, 1997, Vol. 95, 1677-1682.

⁴LCL = 90% exact one-sided lower confidence limit

Table B-3 LIFEPAK CR Plus and LIFEPAK EXPRESS Defibrillator SAS Performance Table for Pediatric ECGs (Continued)

Rhythm Class	ECG Test ¹ Sample Size	Performance Goal ²	Observed Performance Sensitivity or Specificity [LCL] ³
Nonshockable: all other rhythms	433	>95% specificity	99.3% [98.5%]
Intermediate: fine VF	4	Report only	100.0% [56.2%] sensitivity
Intermediate: other VT	7	Report only	42.9% [17.0%] specificity

¹ From Physio-Control ECG database.

Control of Shock Therapy

The Shock Advisory System causes the defibrillator to charge automatically when it detects the presence of a shockable rhythm. When a shockable rhythm is detected, the defibrillator automatically delivers a shock or instructs the user to deliver the shock by pressing the shock button.

Motion Detection

The Shock Advisory System detects victim motion independent of ECG analysis. A motion detector is designed into the defibrillator. MOTION DETECTION can be configured to be ON or OFF.

A number of activities can create motion, including CPR, rescuer movement, patient movement, vehicle movement, and some internal pacemakers. If variations in the transthoracic impedance signal exceed a maximum limit, the Shock Advisory System determines that patient motion of some kind is present. If motion is detected, the ECG analysis is inhibited. The operator is advised by a displayed message, a voice prompt, and an audible alert. After 10 seconds, if motion is still present, the motion alert stops and the analysis always proceeds to completion. This limits the delay in therapy in situations where it may not be possible to stop the motion. However, the rescuer should remove the source of motion whenever possible to minimize the chance of artifact in the ECG.

There are two reasons why ECG analysis is inhibited when the motion alert occurs, and why the rescuer should remove the source of the motion whenever possible:

- Such motion may cause artifact in the ECG signal. This artifact may occasionally cause the Shock Advisory System to reach an incorrect decision.
- The motion may be caused by a rescuer's interventions. To reduce the risk of inadvertently shocking
 a rescuer, the motion alert prompts the rescuer to move away from the patient. This will stop the
 motion and ECG analysis will proceed.

² Automatic External Defibrillators for Public Access Defibrillation: Recommendations for Specifying and Reporting Arrhythmia Analysis Algorithm Performance, Incorporating New Waveforms, and Enhancing Safety. AHA Task Force on Automatic External Defibrillation, Subcommittee on AED Safety and Efficacy. Circulation, 1997, Vol. 95, 1677-1682.

³LCL = 90% exact one-sided lower confidence limit.

APPENDIX C USER'S CHECKLIST

This User's Checklist may be reproduced.

LIFEPAK CR® PLUS LIFEPAK EXPRESS® DEFIBRILLATOR USER'S CHECKLIST



Unit Serial Number	
Department/Location	

		Recommended	Date				
	Instruction	Corrective Action	Initials				
1	Check readiness display for:						
	OK indicator	None.					
	CHARGE-PAK indicator	Replace CHARGE-PAK™ Battery Charger and QUIK-PAK™ Electrode Packet.					
	ATTENTION indicator	Refer to operating instructi	ons.				
	WRENCH indicator	Contact authorized service personnel.					
2	Check Use By date on all Electrode Packets.	Replace electrode packet and CHARGE-PAK if date passed.					
3	Check additional supplies.	Replenish as needed.					
4	Check defibrillator for:						
	Damage or cracks	Contact authorized service personnel.					
	Foreign substances	Clean the device.					
5	Other:						

APPENDIX D

ELECTROMAGNETIC COMPATIBILITY GUIDANCE

Electromagnetic Compatibility Guidance

Table D-1 Guidance and Manufacturer's Declaration - Electromagnetic Emissions

Guidance and Manufacturer's Declaration – Electromagnetic Emissions

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are intended for use in the electromagnetic environment specified below. The customer or the user of the defibrillator should ensure that the defibrillator is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The defibrillator uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The defibrillator is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies
Harmonic emissions IEC 61000-3-2	Not Applicable	buildings used for domestic purposes.
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not Applicable	

Essential Performance

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators maintain safe and effective performance of the defibrillation therapy and patient monitoring functions when operated in the electromagnetic environment specified in Tables 2 through 4.

Limitations Affecting Immunity to Electromagnetic Disturbances

The level of protection from electromagnetic disturbances is limited by several factors, including requirements for protection from third-party defibrillators, patient safety isolation, and maintenance of adequate signal-to-noise ratios for processing patient signals.

Electromagnetic Compatibility Guidance

Table D-2 Guidance and Manufacturer's Declaration - Electromagnetic Immunity

Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are intended for use in the electromagnetic environment specified below. The customer or the user of the defibrillator should ensure that the defibrillator is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±8 kV contact ±15 kV air	The defibrillator is suitable for use in a dry environment.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	Not Applicable	Not Applicable
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	Not Applicable	Not Applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	$<5\% U_{\rm T}$ $(>95\% \ dip \ in \ U_{\rm T})$ for 0.5 cycle $40\% \ U_{\rm T}$ $(60\% \ dip \ in \ U_{\rm T})$ for 5 cycles $70\% \ U_{\rm T}$ $(30\% \ dip \ in \ U_{\rm T})$ for 25 cycles $<5\% \ U_{\rm T}$ $(>95\% \ dip \ in \ U_{\rm T})$ for 5 s	Not Applicable	Not Applicable
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Note: U_T is the a.c. mains voltage prior to application of the test level.

Electromagnetic Compatibility Guidance

Table D-3 Guidance and Manufacturer's Declaration - Electromagnetic Immunity

Guidance and Manufacturer's Declaration – Electromagnetic Immunity				
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance	
environment spe	The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are intended for use in the electromagnetic environment specified below. The customer or the user of the defibrillator should ensure that the defibrillator is used in such an environment.			
			Portable and mobile RF communications equipment should be used no closer to any part of the defibrillator, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.	
			Recommended separation distance	
Conducted RF	3 Vrms	Not Applicable	Not Applicable	
IEC 61000-4-6	150 kHz to 80 MHz outside ISM bands*			
	10 Vrms	Not Applicable	Not Applicable	
	150 kHz to 80 MHz in ISM bands*			
Radiated RF	10 V/m	10 V/m	$d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz	
IEC 61000-4-3	80 MHz to 2.5 GHz	80 MHz to 870 MHz, 900 MHz to 2.5 GHz	$d = 2.3 \sqrt{P}$ for specified frequencies in the range 800 MHz to 2.5 GHz	
		3 V/m	$d = 7.7 \sqrt{P}$ 870 MHz to 900 MHz	
		870 MHz to 900 MHz		
			Where P is the maximum output power rating of the transmitter in watts (w) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).	
			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,‡ should be less than the compliance level in each frequency range.	
			Interference may occur in the vicinity of equipment marked with the following symbol:	
			((🕒))	

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz.

[†]The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 is used in calculating the recommended separation distance for transmitters in these frequency ranges.

Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitter, an electromagnetic site survey should be considered. If the measured field strength in the location in which the defibrillator is used exceeds the applicable RF compliance level above, the defibrillator should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the defibrillator.

Table D-4 Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators

Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators

The LIFEPAK CR Plus and LIFEPAK EXPRESS defibrillators are intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the defibrillator can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the defibrillator as recommended below, according to the maximum output power of the communications equipment.

	Separation distance according to frequency of transmitter m			
Rated maximum output power of transmitter	150 kHz to 80 MHz	80 MHz to 800 MHz	80 MHz to 870 MHz, 900 MHz to 2.5 GHz	870 MHz to 900 GHz
W		d = 1.2√P	d = 2.3√P	d = 7.7√P
0.01	Not Applicable	0.12	0.23	0.77
0.1	Not Applicable	0.38	0.73	2.43
1	Not Applicable	1.2	2.3	7.7
10	Not Applicable	3.8	7.3	24.3
100	Not Applicable	12	23	77

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz.

Note 3: An additional factor of 10/3 is used in calculating the recommended separation distance for transmitters in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas.

Note 4: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

INDEX

Α

About Automatic External	Deleting Data 4-2	Fibrillation 1-3
Defibrillators 1-2	Previous Record 4-2	
About the Defibrillators 1-9	Storing Data 4-2	G
Accessories 5-8	Defibrillation 1-2, 1-3	Getting Started 2-1
ADAPTIV 1-10	Defibrillator	· ·
AED 1-2, 1-3	Operator Settings	Н
Attention Indicator 2-4, 5-2	Motion Detection 6-4	Heart Attack 1-3
В	Device Identification 6-2	
Basic Steps 3-3	Disconnecting Electrode	1
Blue Plastic 2-6, 2-7	Pads 3-5	Impedance 1-3
Diue Plastic 2-6, 2-7		Indications for Use 1-2
С	E	Indicators 2-4
	ECG 1-3	Inside Features 2-5, 2-6
Cardiac Arrest 1-3	Electrode Connector 2-6, 2-7	Inspecting the Defibrillator 2-4
Caring for the Defibrillator 5-1	Electrode Indicators 2-6, 2-7	IrDA Port 1-10, 2-4, 2-5
Carrying Handle 2-4, 2-5	Electrode Packet 2-5, 2-6	_
CHARGE-PAK 1-10, 2-4, 2-5	Electrode Packet Anchor Pin 2-5,	J
CHARGE-PAK Indicator 2-4, 5-2	2-6	Joule 1-3
Cleaning Methods 5-3	Electrode Packet Release	
Cleaning the Defibrillator 5-2	Handle 2-5, 2-6	L
Clinical Summary A-5	Electrode Pad	Labels 2-4
Continuous ECG 4-2	Recycling/Disposal 5-8	LED 1-3
Controls 2-4	Electrode Pads 2-6, 2-7	Lid 2-4
CPR 1-3	Energy Protocol 6-2	LIFEPAK Defibrillator
CPR Time 6-3	Energy Sequence	Accessories 1-10
D	Energy Levels 6-2	Automated Operation
D	Environmental Specifications A-3	Fully Automatic 1-10
Data Storage 4-2	Event and Test Log 4-3	Semiautomatic 1-10
Data Storage Characteristics A-4	Event and Test Log Reports 4-3	Automatic Self-Test 1-10
Data Stored by the	Event Log 4-2	Capabilities and
Defibrillator 4-2	ŭ	Features 1-10
Data Transfer		realures 1-10

Current Record 4-2

F

Customized Setup 1-10	Q	Voice Prompts and Tones 3-5
Data Management 1-10	Quick Reference Card 2-5, 2-6	Voice Prompts,
Defibrillation	QUIK-PAK Electrodes 1-10	Troubleshooting 3-6
Electrodes 1-10	R	W
Defibrillation Waveform 1-10	Readiness Display 2-4	Warnings and Cautions 1-5, 3-2
Heart Rhythm Analysis 1-10	Readiness Display,	General 1-5
Operator Settings Motion Detection 1-11	Troubleshooting 3-6	Warranty Information 5-8
Power System 1-11	Recycling Information 5-8	What to Do After Emergency
Readiness Display 1-11	Regular Maintenance 5-2	Medical Personnel
Troubleshooting 3-6	Reports Stored 4-2 Responder 1-3	Arrive 3-5 What to Do After Using the
Log Reports 4-3	Responding to a Sudden Cardiac	Defibrillator 3-5
5 .	Arrest Victim 3-2	Where to locate your
M		Defibrillator 2-3
Maintaining a State of	S	Why the Need for
Readiness 5-2	Safety Information 1-4	Defibrillators 1-2
Managing Data 4-1 Motion Detected 3-6	Safety Terms 1-4 Safety Warnings 2-4, 2-5	Wrench Indicator 2-4, 5-2
Motion Detection 6-4	SAS 1-3	
Myocardial Infarction 1-3	Serial Number Label 2-4, 2-5	
	Service 5-8	
N	Shock Advisory System B-1	
Nonshockable Rhythm 1-3	Shock Button 2-6, 2-7 Shockable Rhythm 1-3	
0	Speaker 2-6, 2-7	
Obtaining Authorized Service 5-8	Specifications A-1	
OK Indicator 2-4, 5-2	Stack Shocks 6-3	
ON-OFF Button 2-4	Storing the Defibrillator 5-2	
Opening the Electrode	Summary 4-2	
Packet 3-3	Supplies 5-8 Symbols 1-7	
Operating Settings 6-2 About 6-1	Symbols 1-7	
Operating Settings and Setup	Т	
Configuration 6-2	Terminology 1-3	
Operator Settings	Test and Service Data 4-2	
CPR Time Settings 6-3	Text Conventions 1-4	
Device Date 6-3	Time Zone 6-4 Training Tools 5-8	
Device ID 6-2	Treating a Victim 3-3	
Device Time 6-3	Troubleshooting 3-6	
Energy Protocol 6-2	Troubleshooting During Victim	
Energy Sequence 6-2 Pulse Check 6-3	Use 3-6	
Pulse Prompt 6-3	U	
Stack Shocks 6-3	Unpacking and Inspecting 2-2	
Time Zone 6-4	Use By Date 2-5, 2-6	
Voice Prompt Volume 6-3	User 1-3	
Outside Controls Indicators and	User Interface Specifications A-2	
Labels 2-4	Using the Defibrillator 3-1, 3-3	
Overview of Data Storage 4-2	V	
P	V Ventricular Fibrillation 1-2, 1-3	
Physical Characteristics A-3	Ventricular Tachycardia 1-3	
Placing Electrode Pads 3-3, 3-4,	Victim 1-3	
3-6	Victim Care	
Problem Conditions 3-6	Transferring 3-5	
Pulse Check 6-3	Treating the Victim 3-2	
Pulse Prompt 6-3	Troubleshooting 3-6	
	Voice Prompt Volume 6-3	













